

Mojia


CHARIS

Microcurrent

User Manual



Getting To Know Your Device

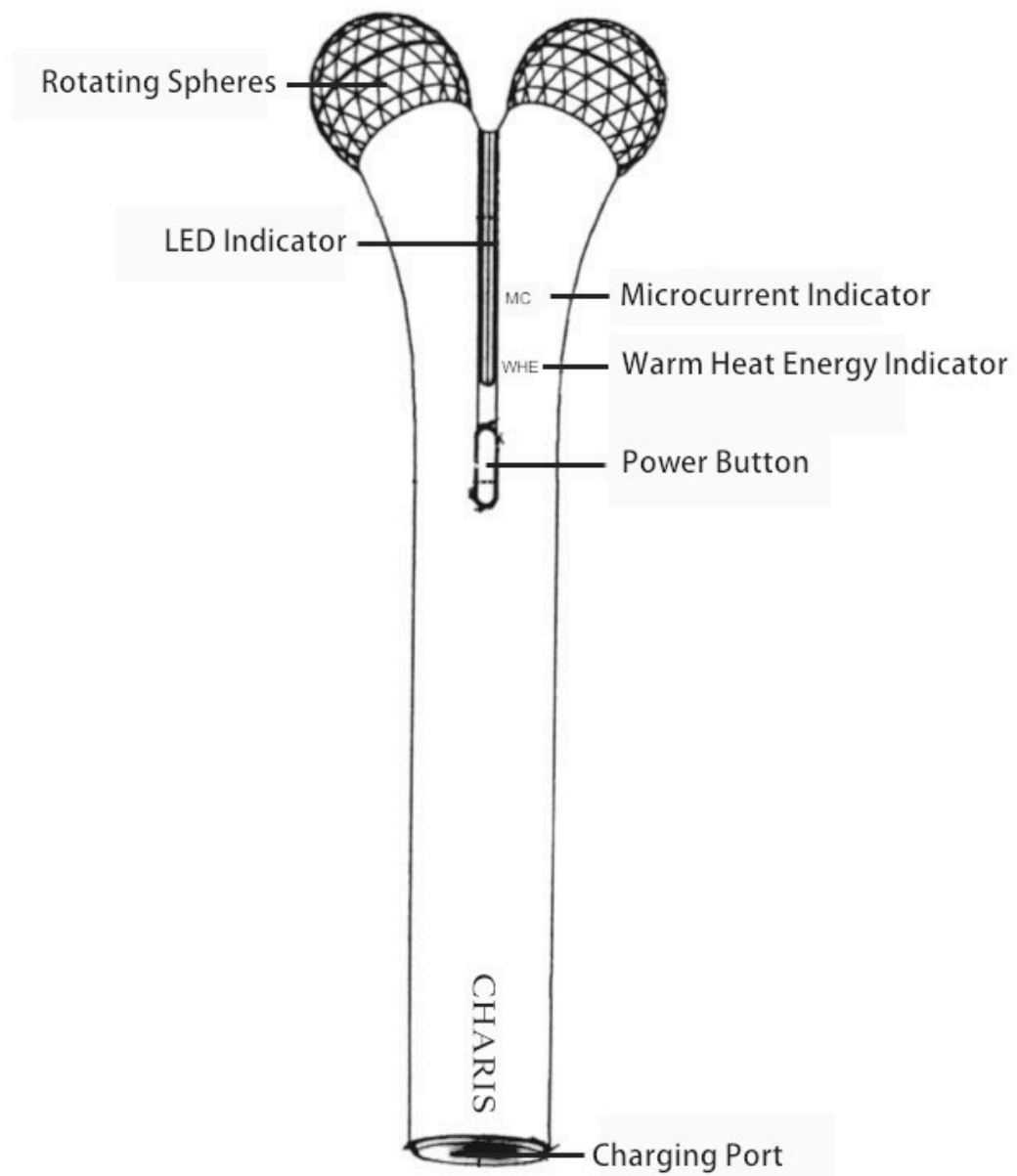

Before Your First Use

- When you receive your Mojia CHARIS you will need to charge the battery for 3 hours before turning it on.
- Do not use the device while charging.
- The LED Indicator strip will pulsate during charging.
- When the battery is fully charge the LED Indicator strip will remain solid.
- Press and hold the power button for 3 seconds to turn the device on.
- By default the WHE (warm heat energy) light comes on.
- Press the power button once for the MC mode.
- Press the power button to increase the MC intensity. This device has three settings.
- Press and hold the power button for 3 seconds to turn the device off.

The Mojia CHARIS can be used in the WHE mode or the WHE and MC mode.

- The device switches off automatically after 6 minutes.
- The LED Indicator strip will flash 3 times when battery is low.

The Mojia CHARIS can be used as a massage tool, using the WHE mode only you can glide or use circular movements for relaxing muscle tension.

It can also be used to increase product penetration. This mode takes approximately 20 seconds to start feeling warm.

If your skin needs a quick pick-me-up, de-puffing, or soothing, place the Mojia CHARIS in the fridge for a maximum of 10 minutes and treat your skin to some cool relief. Do NOT turn the device on.

The MC mode incorporates WHE (warm heat energy). This combination improves skin and muscles response to the micro-currents as well as product absorption to the skin.

Before Any Treatment

It is important to make sure your skin is oil free and clean.

- Wash face and neck area with an oil free cleanser and pat dry.
- Make sure all makeup and residue is completely removed.

Important- both rotating spheres must be in contact with the skin surface.



STEP 1

Apply conductive gel liberally to one side of the neck and jawline. Place the spheres at the base of the neck. Make sure both spheres are touching the skin. Slowly glide upward to the jawline. Once reached hold for 3 seconds. Repeat 3 times. Move to the next section and repeat.

Always avoid the thyroid area, as shown in red.



STEP 2

Follow the diagram, moving outwards.

Hold for 3 seconds. Starting at the base of the neck. Repeat 3 times. Move up to the next section and repeat.

Continue the same movements starting in the middle under the chin and jawline, moving outwards.

Repeat 3 times.



STEP 3

Apply conductive gel to one side of the jaw and cheek area. Start with one sphere under the jawline and the other on top, so the jawline is wedged between the spheres. Slowly glide outward. Hold for 3 seconds. Repeat 3 times.


Move up to the next area and repeat.



STEP 4

Starting at the jawline near the corner of your mouth slowly glide up to the cheek bone. Hold for 3 seconds. Repeat 3 times. Move to the next area and repeat.

Repeat Steps 3 & 4 to the opposite side.

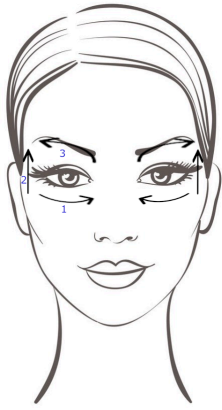

STEP 5

Gently apply conductive gel around the eye area as indicated in the diagram. From gently and slowly glide from the outer corner of the eye toward the nose. Do not hold. Repeat 3 times. Glide upwards over the temple area (crows feet) as indicated by the arrow. Hold for 3 seconds and repeat 3 times.


With one sphere below the brow and one above, slowly glide outwards toward the hairline. Hold for 3 seconds and repeat 3 times.



STEP 6


Placing one sphere at the end of the brow and the other near the middle, gently lift and hold for 3 seconds. Release the lift, then gently lift again and repeat 3 times.

Repeat Steps 5 & 6 to the opposite side.

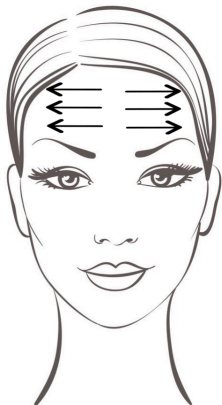

STEP 7

Apply the conductive gel to one side of the forehead. Gliding slowly outwards from the middle of the forehead towards the hairline. Hold for 3 seconds. Move to the next area and repeat.

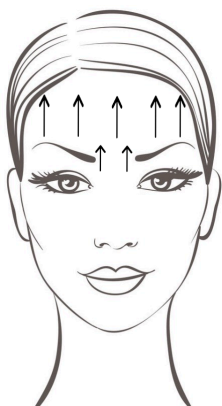

STEP 8

Starting at the top of the brow slowly glide upwards. Hold for 3 seconds then move to the next area.

For the 11' s (in between the brows) roll the spheres slowly over the area for approximately 10 seconds.

Repeat Steps 7 & 8 on the opposite side.

Conductive Gel

Mojia Australia recommends using our Mojia Electra Conductive Booster Gel which can be left on and has been specially formulated with 17 active ingredients for maximum results.



Cleaning your device

Always clean the Mojia CHARIS after every use.

Do not immerse in water.

Clean the rotating spheres with a damp cloth and wipe dry.
Store device away from any heat source, direct sunlight or moisture.

Precautions and Warnings

- Do not use for any other purpose than its intended use as recommended in this manual.
- Do not use if the device is defective as this can lead to personal injury.
- Do not use directly on the eyes.
- Do not use on infected, inflamed, damaged or broken skin.
- Stop using the device if you experience an abnormal reaction.
- A slight tingling sensation may be felt, this is normal and no cause for concern. You can decrease the level of intensity to reduce this sensation.
- Do not use if you suffer from Seizures or Epilepsy.
- Do not use if you suffer from heart issues or have a pacemaker.
- Wait two weeks after having botox or fillers before using this device.
- Do not use the device if you have an allergy to metals.
- Please consult your medical practitioner if unsure or have any concerns.
- Keep out of reach of children.
- Do not use the device whilst charging.
- Do not modify, disassemble and/or repair the device. Please contact Mojia for assistance.

Warranty

1-YEAR LIMITED WARRANTY

Mojia warrants this device for a period of ONE (1) year after the date of purchase. The warranty covers defects due to faulty workmanship. It covers working parts that affect the function to the device. It does not cover cosmetic deterioration caused by normal wear and tear, damaged caused by accident, misuse or neglect.

Please refer to the Terms & Conditions of Mojia Australia.

Disclaimer: Users of this device do so at their own risk. Mojia Australia do not assume any responsibility or liability for any injuries or damages resulting from the use of this device, whether it be physical or otherwise.



European
Conformity



Restriction Of
Hazardous
Substances



Do not dispose in
household bin. Take to
recycling point.

Mojia

CHARIS
Microcurrent

Mojia Australia

www.mojiaaustralia.com

info@mojiaaustralia.com