


✦beautiva™

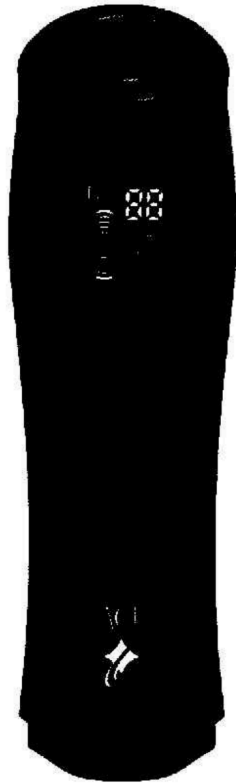

User Manual

Be sure to read the user manual before use.

Precautions & Warnings

Please read the user manual carefully before using the device for the first time.

SAFETY

- Only use this device for its intended purpose.
- Only use this device with the power cord provided.
- Plug into power outlet and into device before switching power outlet on.
- Always turn the device off and then plug from the power outlet before cleaning.
- Please stop use immediately if any skin abnormalities occur and seek medical advice.
- Do not use this device if it is cracked or damaged.
- Do not use if the power cord is damaged.
- Do not touch the device should it fall into water.
- Do not place or store the device near any source of water, such as a sink.
- Do not submerge the device under water.
- When plugged in, do not leave the device unattended.
- Keep the device out of the reach of children.
- Do not use the device over eyelids, thyroid area, ears, breasts or genital areas.

DO NOT USE IF YOU HAVE ANY OF THE FOLLOWING

- Heart disease, abnormal blood pressure, malignant tumours, cancer, skin cancer.
- Migraines.
- Seizures.
- Pacemaker.
- Have any active skin conditions, such as, but not limited to, dermatitis, psoriasis, eczema, rash or open sores.
- Have had recent chemical peels, laser resurfacing or similar treatments.
- Botox or filler within 2 weeks.
- Sunburn.

If unsure please contact your practitioner for further advice.


CLEANING AND MAINTENANCE

- Turn the device off and unplug from the power outlet.
- Use a damp cloth to clean with a damp cloth making sure to remove any product in between the probes.
- Dry thoroughly.

PRODUCT STORAGE


- Keep the device unplugged when not in use.
- Do not place in direct sunlight.
- Store device in its box in a dry place.
- A certain amount of heat when using the device is normal, however, if you feel extreme or abnormal heat cease use immediately.
- Do not disassemble, modify or repair this device.
- Do not submerge device in water or get it wet, this may lead to serious injury or electric shock. It may also cause the device to malfunction.

Product & Accessories

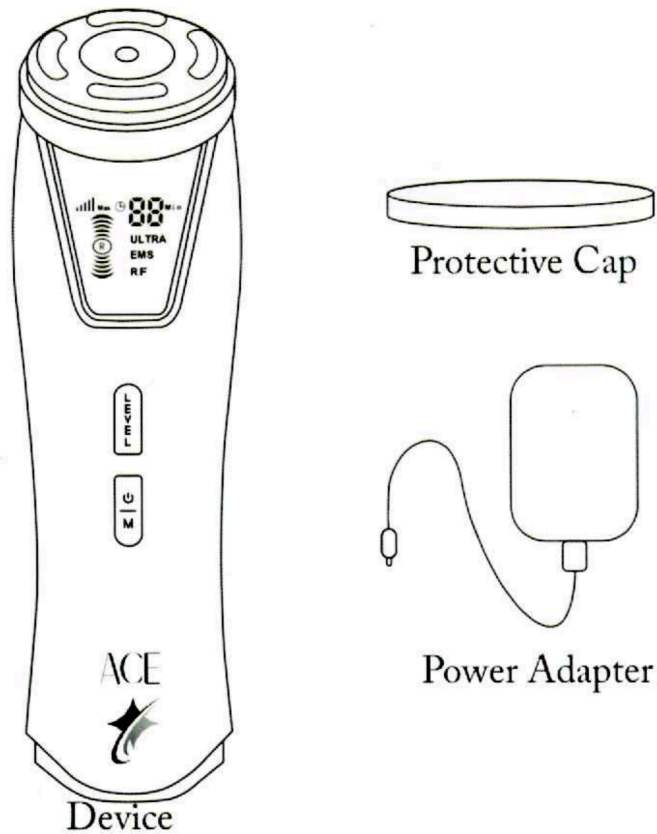

Power Adapter
Enter: 100V-240V~50/60Hz 0.6A Output: 24V 1A

Get To Know Your ACE

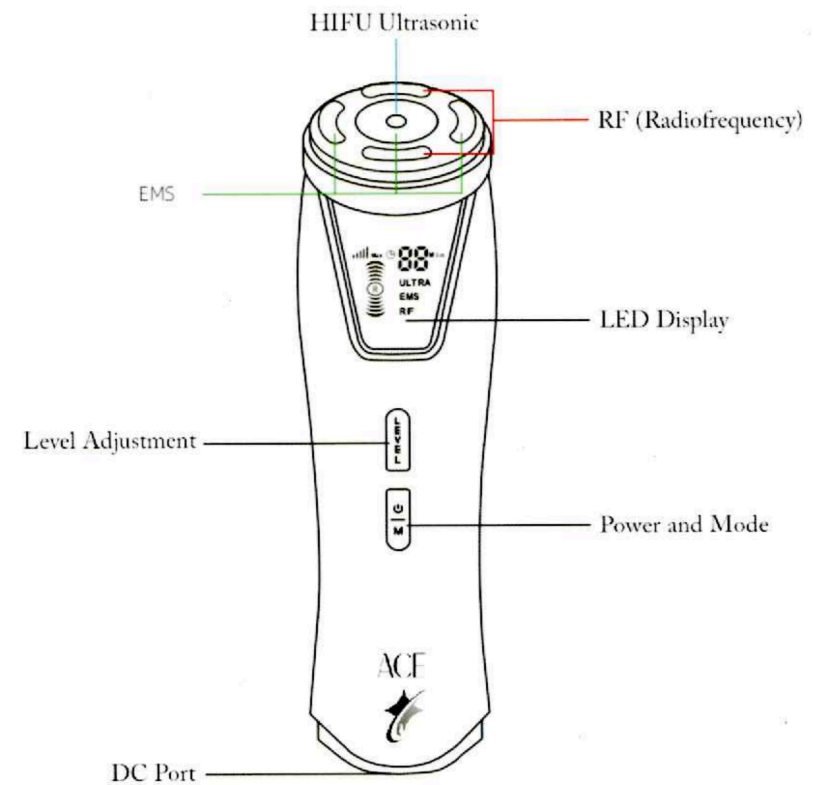


Product Description

The Beautiva Ace is a revolutionary device that offers a combination of hifu, radiofrequency, and EMS treatments all in one device. What makes this device unique, is its ability to use radiofrequency and EMS simultaneously, providing users with an efficient, effective treatment.


This advanced technology also allows for a faster treatment time, making it the perfect choice for those looking for a quick yet effective solution.




HIFU (High Intensity Focused Ultrasound) uses ultrasound energy to target the deeper layers of the skin, stimulating collagen production and tightening the skin for a more youthful appearance. The energy is delivered to the skin without damaging the surface.

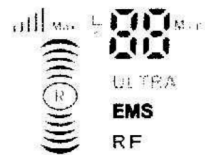

EMS (Electrical Muscle Stimulation) uses mild current to stimulate the muscle by sending electrical impulses that cause it to contract.

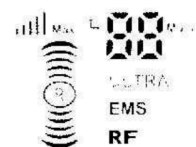

RF (Radiofrequency) uses low radiofrequency waves to deliver heat deep into the dermis to stimulate fibroblasts to create collagen and elastin.

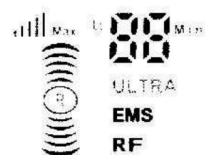

EMS/RF - this type of treatment is known as Combination Therapy. It involves the use of both radiofrequency and EMS (electrical muscle stimulation) to target deep layers of the skin and muscle tissue, resulting in a more comprehensive and effective treatment. The radiofrequency energy heats the skin and underlying tissue, stimulating collagen production and improving skin elasticity and texture, while the EMS helps to tone and lift the skin and reduce wrinkles.

How to Use Your Device

1. Plug in the power supply.
2. Press the power button until the device turns on. The LED screen will light up. The default mode is ULTRA (HIFU Ultrasound). The default intensity level is 1. There are 5 intensity levels for each mode. The timer is set at 10 minutes.
3. Toggle between modes with a quick short press of the mode button.
4. A beep will sound with each press.
5. Toggle between intensity levels by pressing the button. 1 being the lowest and 5 being the highest.
6. BLUE LED light indicates HIFU mode.
7. GREEN LED light indicates EMS mode.
8. RED LED light indicates RF mode.
9. Alternating GREEN and RED LED lights indicates combined EMS and RF mode.
10. When using combined EMS and RF mode increasing the intensity levels only increased the RF. The EMS intensity is set and cannot be increased.
11. To turn the device off, press and hold the power button until device turns off.

NOTE

- A beep will sound each time a button is pressed.
- The 10 minute timer begins to countdown when the device is turned on.

Treatment


Preparation

- Use an oil free cleaner to remove any makeup, oils or dirt. Pat dry.
- Remove any jewellery that may come in contact or interfere with the device.
- During your first treatment it is recommended to start on the lowest intensity level and adjust to your skins tolerance.

HIFU

The Ace is designed with advanced technology and cooling plates so it doesn't get hot in HIFU mode.

- Divide your face into 5 treatment sections
 - Left cheek
 - Right cheek
 - Forehead
 - Under chin
 - Neck

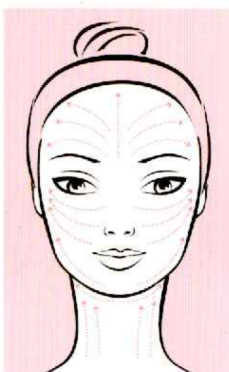

Apply a liberal amount of gel to the section you are working on at a time.

HIFU - Glide very slowly in the direction of the arrows. Each section should take approximately one minute when doing one pass (each arrow), however, you can do two passes on problem areas such as jawline and under the chin.

- * A slight tingle sensation may be felt.
- * Never use over the thyroid or on the eyelids.

DO NOT EXCEED A 10 MINUTE TREATMENT
Treatment - ONCE A WEEK for 10 minutes.

EMS

- Divide your face into 3 treatment sections.
 - Left cheek
 - Right cheek
 - Under Chin



Apply a liberal amount of gel to the section you are working on at a time.

Use gliding movements in the direction on the arrows.

- * Never use over the thyroid or on the eyelids.

DO NOT EXCEED A 10 MINUTE TREATMENT

Treatment - Up to 5 times weekly for 10 minutes.

RF

- Divide your face into 5 treatment sections.
 - Left cheek - lower and upper
 - Right cheek - lower and upper
 - Under Chin



Apply a liberal amount of gel to the section you are working on at a time.

Using continuous small circular movements on each area for 2 minutes.


Treatment on the forehead is recommended at a low intensity.

Do not hold in one place.

- * Never use over the thyroid or on the eyelids.

DO NOT EXCEED A 10 MINUTE TREATMENT

Treatment - Up to 3 times weekly for 10 minutes.

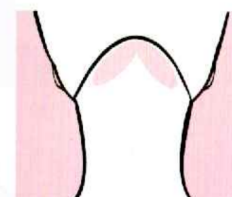


RF & EMS Combined therapy

- Divide your face into 5 treatment sections.
- Left cheek - lower and upper
- Right cheek - lower and upper
- Under Chin

Using either circular or gliding movements

* Never use over the thyroid or on the eyelids.

RF & EMS Combined therapy provide quick 5 minute treatments that can be used daily.

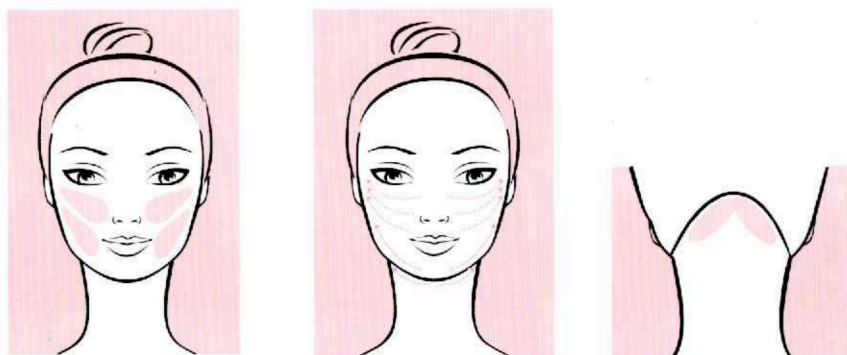

Important:

- Do not use Hifu and RF on the same day.
- To protect the device from damage it will turn off if overused.
- The device may feel warm when holding and this is normal.
- Skin may appear pink after treatment and this will subside quickly.

Technical Information

Functions	HIFU, RF, EMS and LED
HIFU	5J/cm ²
Work area	4cm ²
Frequency	3MHz
Depth	Up to 4.5mm
RF Power	4J/cm ²
Frequency	2MHz
Depth	3mm
Power	36W
Power source	15V AC adapter
Voltage	110/240V, 50/60 Hz biovolt
Weight	145gm



